

Various magnifications (from 2.3x to 6x), working distances, wearing modes (unique S-FRAME, Lightweight Headband and Headband Professional L) and LED LoupeLight illumination systems are available. They provide homogeneous, reflex-free illumination even in deep cavities. State-of-the-art LED HQ technology and design!

HEINE Binocular Loupes are used in many disciplines: e.g. General Medicine, Dermatology, Surgery, Cosmetic Surgery, Neurosurgery, ENT, Ophthalmology, Veterinary Medicine, Dentistry etc.

Definitions:

∴ Working distance (x)


The distance, at which the optics of a loupe are sharply-focused.

∴ Depth of field (y)

The range over which the image remains sharply-focused.

∴ Field of view (z)

The area that can be seen under magnification at the nominal working distance.

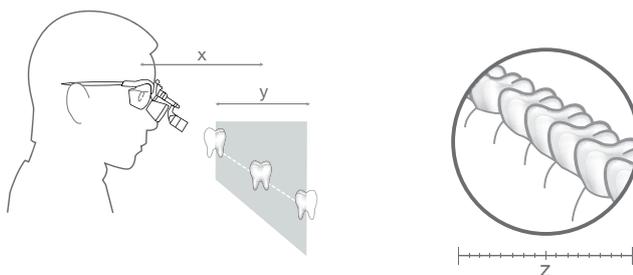

How to choose the correct Loupe:

∴ Choose your magnification

Choose the lowest magnification that meets your needs. The lower the magnification, the larger the field of view and the easier the loupe will be to use.

∴ Choose your optics

HR: High-Resolution achromatic optics with excellent optical qualities.

HRP: High Resolution Prismatic optics with excellent optical qualities for 3.5x magnification and higher.

C: Economy.

∴ Choose your working distance

This is an individual choice depending on the examiner's size and working position. The longer the working distance, the bigger the field of view.

∴ Choose your mounting option

Spectacle frame or Headband.

∴ Add illumination

Coaxial bright light is especially-important when using magnification.

HEINE Binocular Loupes

	HR	HRP	C
			C 2.3K C 2.3
Magnification			
2.3x			<input type="checkbox"/>
2.5x	<input type="checkbox"/>		
3.5x		<input type="checkbox"/>	
4x		<input type="checkbox"/>	
6x		<input type="checkbox"/>	
Working distance			
250mm (with close-up lenses)	<input type="checkbox"/>	--	--
340mm	<input type="checkbox"/>	<input type="checkbox"/>	<input type="checkbox"/>
420mm	<input type="checkbox"/>	<input type="checkbox"/>	--
450mm	--	--	<input type="checkbox"/>
520mm	<input type="checkbox"/>	--	--
Mounting options			
S-FRAME			
	<input type="checkbox"/>	<input type="checkbox"/>	--
Lightweight Headband			
	--	--	<input type="checkbox"/>
Professional Headband			
	<input type="checkbox"/>	<input type="checkbox"/>	--
Simple Metal Frame			
	--	--	<input type="checkbox"/>
Accessories			
Close-up Lenses*	<input type="checkbox"/> (for 340mm)	--	--
Protective Lenses	<input type="checkbox"/>	<input type="checkbox"/>	<input type="checkbox"/>
LED LoupeLight (Pages 113/114)			
	<input type="checkbox"/>	<input type="checkbox"/>	--
Features			
i-View Function**	<input type="checkbox"/>	<input type="checkbox"/>	--
Modular	<input type="checkbox"/>	<input type="checkbox"/>	--
Water-resistant	<input type="checkbox"/>	<input type="checkbox"/>	--
Achromatic Optics	<input type="checkbox"/>	<input type="checkbox"/>	--
Optional Prescription Correction	<input type="checkbox"/>	<input type="checkbox"/>	--
Page	106 – 107	108 – 109	110



* To reduce the working distance to 250mm.

** **i-View function:** Provides any angle of view and flips up the optics independently of the illumination beam.



HEINE® HR 2.5x High Resolution Binocular Loupes with i-View loupe mount



∴ Multi-coated Achromatic Optics

∴ Exceptional Focal Range

∴ Crisp, High-Resolution Image

∴ Unique i-View Loupe Mount

HR 2.5x Optics with i-View

- ∴ **Achromatic optics** for crisp, bright images with colour-correction.
- ∴ **Lightweight** (51 g).
- ∴ **Water-resistant.**
- ∴ **Independent adjustment** of left and right optics for perfect PD adjustment.
- ∴ **High magnification** with a large field of view and excellent depth of field.
- ∴ **i-View loupe mount** provides any angle of view and flips up the optics independently of the optional LED LoupeLight.

HR 2.5x Specifications

Working distance in mm	340	420	520
Field of view in mm	90	130	180
Depth of field in mm	110	180	260



Optional clip-in correction frame

Mounting Options:

∴ HEINE S-FRAME®

- ∴ **Maximum comfort** with a modern, ergonomic design.
- ∴ **Malleable temples and nose pad** can be customized to fit any face perfectly.
- ∴ **Innovative materials** – Carbon-Kinetium reduces weight (24 g).
- ∴ **Tough polycarbonate protective lenses.**
- ∴ **Optional clip-in correction frame** for prescription lenses.



Professional L Headband with S-GUARD

∴ HEINE® Professional L Headband **with** S-GUARD®

- ∴ **S-GUARD splash protection.**
- ∴ **Ergonomic nose pad** for comfort during long procedures.
- ∴ **Fully-adjustable for height and width to fit any head.**
- ∴ **Hygienic and comfortable leather cushions.**
- ∴ **Optional clip-in correction frame** for prescription lenses.



Professional L Headband without S-GUARD

∴ HEINE® Professional L Headband **without** S-GUARD®

- ∴ **Ideal for spectacle-wearers in combination with their own spectacles.**
- ∴ **Ergonomic nose pad** for comfort during long procedures.
- ∴ **Fully-adjustable for height and width to fit any head.**
- ∴ **Hygienic and comfortable leather cushions.**

HEINE® HR 2.5x High Resolution Binocular Loupes
with i-View loupe mount



Set A



Set B
with S-GUARD



Set C
without
S-GUARD

HR 2.5x Sets with i-View loupe mount

Working distance in mm	340	420	520
Set A: HR 2.5x optics with i-View on S-FRAME Protective lenses 2 sterilisable swivel levers Retaining cord Cleaning fluid Cleaning cloth Carrying case	C-000.32.355	C-000.32.356	C-000.32.567
Set B: HR 2.5x optics with i-View on Professional L Headband S-GUARD splash protection Protective lenses 2 sterilisable swivel levers Cleaning fluid Cleaning cloth	C-000.32.365	C-000.32.366	C-000.32.367
Set C (without S-GUARD): HR 2.5x optics with i-View on Professional L Headband Protective lenses 2 sterilisable swivel levers	C-000.32.865	C-000.32.866	C-000.32.867

For details of accessories and LED LoupeLight see pages 112 – 114.



HEINE® HRP High Resolution Prismatic Binocular Loupes with i-View loupe mount



- ⌘ Multi-coated Achromatic Optics
- ⌘ Prismatic for highest magnification
- ⌘ Crisp, High-Resolution Image
- ⌘ Unique i-View Loupe Mount

HRP Optics with i-View

- ⌘ **Achromatic optics** for crisp, bright images from edge-to-edge.
- ⌘ **Lightweight** (85g).
- ⌘ **Water-resistant.**
- ⌘ **Independent adjustment** of left and right optics for perfect PD adjustment.
- ⌘ **Highest magnification:** Choice of 3.5x, 4x and 6x.
- ⌘ **i-View loupe mount** provides any angle of view and flips up the optics independently of the optional LED LoupeLight.

HRP Specifications

Working distance in mm	HRP 3.5x/420	HRP 4x/340	HRP 6x/340
Field of view in mm	65	50	40
Depth of field in mm	60	40	30

Mounting Options:



⌘ HEINE S-FRAME®

- ⌘ **Maximum comfort** with a modern, ergonomic design.
- ⌘ **Malleable temples and nose pad** can be customized to fit any face perfectly.
- ⌘ **Innovative materials** – Carbon-Kinetium reduces weight (24g).
- ⌘ **Tough polycarbonate protective lenses.**
- ⌘ **Optional clip-in correction frame** for prescription lenses.



⌘ HEINE® Professional L Headband **with** S-GUARD®

- ⌘ **S-GUARD splash protection.**
- ⌘ **Ergonomic nose pad** for comfort during long procedures.
- ⌘ **Fully-adjustable for height and width to fit any head.**
- ⌘ **Hygienic and comfortable leather cushions.**
- ⌘ **Optional clip-in correction frame** for prescription lenses.



HEINE® Professional L Headband **without** S-GUARD®

- ⌘ **Ideal for spectacle-wearers in combination with their own spectacles.**
- ⌘ **Ergonomic nose pad** for comfort during long procedures.
- ⌘ **Fully-adjustable for height and width to fit any head.**
- ⌘ **Hygienic and comfortable leather cushions.**

HEINE® HRP High Resolution Prismatic Binocular Loupes
with i-View loupe mount



Set A



Set B
with S-GUARD



Set C
without
S-GUARD

HRP Sets with i-View

Working distance in mm	HRP 3.5x/420	HRP 4x/340	HRP 6x/340
Set A: HRP optics with i-View on S-FRAME Protective lenses 2 sterilisable swivel levers Retaining cord Cleaning fluid Cleaning cloth Carrying case	C-000.32.430	C-000.32.431	C-000.32.432
Set B: HRP optics with i-View on Professional L Headband S-GUARD splash protection Protective lenses 2 sterilisable swivel levers Cleaning fluid Cleaning cloth	C-000.32.440	C-000.32.441	C-000.32.442
Set C (without S-GUARD): HRP optics with i-View on Professional L Headband Protective lenses 2 sterilisable swivel levers	C-000.32.840	C-000.32.841	C-000.32.842

For details of accessories and LED LoupelLight see pages 112 – 114.



HEINE® C 2.3 Binocular Loupes

The cost-effective all-purpose loupe

Ideal for dental laboratories and general medical examination.



C 2.3

- ∴ **Crisp, clear images.** Multi-coated, scratch-resistant silica glass lenses produce crisp, clear images.
- ∴ **Flexibility for precise adjustment.** Independently adjustable PD controls ensure precise alignment of the optics.
- ∴ **Super lightweight and comfortable – only 46g.** Flexible coil temples wrap-around to fit the user’s anatomy securely and comfortably. Spring-loaded temples will adjust to the user’s anatomy. Lightweight metal frame.
- ∴ **The ideal all-purpose loupe.** The feather-light weight, low-profile design and low cost of the C 2.3 loupe make it the ideal all-purpose loupe.
- ∴ **Working distance.** 340mm or 450mm (13 or 18 inches).
- ∴ Also available as C2.3K on Lightweight Headband.



C 2.3K

Specifications	C2.3x/340mm	C2.3x/450mm
Magnification	2.3x	2.3x
Working distance	340mm	450mm
Field of view dia.	110mm	130mm
Weight	46g	46g
Depth of Field	120mm	190mm

User's refractive error	Effective working distance C2.3x/340mm	Effective working distance C2.3x/450mm
+2	440mm	670mm
+1	390mm	560mm
0	340mm	450mm
-1	300mm	420mm
-2	270mm	370mm

Important information for those who wear glasses:

The C 2.3K Series loupe can also be worn without the use of spectacles. Should you prefer to examine without glasses, the charts illustrate how your refractive error will affect the loupe working distance. If you are unsure of your prescription, please consult your eye care provider.

C2.3/C 2.3K Binocular Loupes	
C2.3/340 Loupe in case	C-000.32.039
C2.3/450 Loupe in case	C-000.32.202
C2.3K/340 Headband Loupe in case	C-000.32.230
C2.3K/450 Headband Loupe in case	C-000.32.231

Accessories



C2.3/C2.3K Binocular Loupe Accessories	
Protective case for the C 2.3 Loupes and accessories	C-000.32.038
Detachable protective lenses for C 2.3 and C 2.3K. To protect the optics. Can be sterilized by conventional methods.	C-000.32.097

HEINE Binocular Loupes

Ordering Information



Optics with i-View loupe mount for S-FRAME in carrying case			
Working distance in mm	340mm	420mm	520mm
HR 2.5x	C-000.32.595	C-000.32.596	C-000.32.597
HRP 3.5x HRP 4x HRP 6x	C-000.32.721 C-000.32.722	C-000.32.720	



Optics with i-View and S-GUARD for Professional L Headband			
Working distance in mm	340mm	420mm	520mm
HR 2.5x	C-000.32.585	C-000.32.586	C-000.32.587
HRP 3.5x HRP 4x HRP 6x	C-000.32.711 C-000.32.712	C-000.32.710	



Optics with i-View for Professional L Headband (without S-GUARD)			
Working distance in mm	340mm	420mm	520mm
HR 2.5x	C-000.32.733	C-000.32.734	C-000.32.735
HRP 3.5x HRP 4x HRP 6x	C-000.32.737 C-000.32.738	C-000.32.736	



HEINE Binocular Loupes Accessories

Ordering Information



Mounting options

S-FRAME*

For **HR** and **HRP** Binocular Loupes incl. retaining cord and cleaning fluid

C-000.32.302

* with protective lenses, large (lg).

Professional L Headband

For **HR** and **HRP** optics

C-000.32.549



Accessories for optics

HR Close-up Lenses in case

HR 2.5x/340 mm to 250 mm. The close-up lenses reduce the working distance of the binocular loupes.

Close-up Lenses for **HR 2.5x/340mm** optics

C-000.32.523

Protective Lenses

To protect the optics of the **HR** and the **HRP** loupes

Protective Lenses for **HR** loupes (5 pairs) [01]

C-000.32.537

Protective Lenses for **HRP** loupes [02]

C-000.32.027

Sterilisable Swivel Levers

For **HR** and **HRP** loupes with i-View

C-000.32.551

Carrying Case

For Binocular Loupe optics or for all sets with S-FRAME

C-000.32.552

Accessories for S-FRAME

S-FRAME retaining cord

C-000.32.305

S-FRAME replacement part set

with side rubbers, nose pad (1 piece) and screw

C-000.32.308

Clip-in correction frame

for prescription lenses

C-000.32.309

Protective lenses, large (lg), 1 pair (Standard)

C-000.32.306

Protective lenses, small (sm), 1 pair

C-000.32.307

Cleaning fluid for lenses, pack of 2

C-000.32.542

HEINE® LED LoupeLight

Compact LED LoupeLight – now in LED^{HQ}



- ⚡ Compact and Lightweight
- ⚡ Coaxial Illumination
- ⚡ Up to 40.000 Lux/250mm distance




- ⚡ **LED in HQ** – the new standard in LED illumination defining optimal light intensity, homogeneity and colour rendering for each individual examination situation.
- ⚡ **Compact, lightweight (24 g)**. Ensures perfect comfort.
- ⚡ **Aluminium Housing**. Exclusive aluminum housing protects the optics against dust and moisture to ensure an extraordinary long lifespan.
- ⚡ **Coaxial Illumination**: Compact coaxial design ensures a completely shadow-free image and allows for excellent illumination of difficult to see areas.
- ⚡ **Individual adjustment of the angle of view**. The right angle for every situation.
- ⚡ **Optional Filter**. Polarisation filter P2 improves contrast and yellow filter reduces blue light.


Technical Specification

Illumination	Up to 40,000 Lux at a working distance of 250 mm
Colour Temperature	4,500 Kelvin
LED Operating Life	Up to 50,000 hours
Spot Size	The illuminated field is designed to complement the HR/HRP loupes. 81 mm dia. at 420 mm, 95 mm dia. at 520 mm working distance
Angle of Declination	Stepless and individually adjustable for every examination
Operating Time	17 hours continuous ON time at full power with mPack LL

- ⚡ **Precision multiple lens optical system** consisting of four high-quality lenses. Completely homogeneous illumination of the whole field of view and a clearly defined light spot.

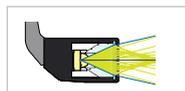

Conventional System [01]

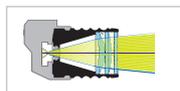

HEINE Precision Lens System [02]

HEINE® mPack LL

Power Source



- ⚡ **Portable power pack**. Does not tie the user to desk, wall or static power supplies.
- ⚡ **Li-ion Technology**. Triple capacity compared with current NiCd systems. No memory effect, autonomy for 17 hours of continuous ON time at full power.
- ⚡ **Fast Charge**. Charges to full capacity in 2 hours.
- ⚡ **Charge Status Display**. Always lets you know how much charge you have left.
- ⚡ **Stepless light intensity control**. Optimum brightness setting prevents reflexes.
- ⚡ **Flexible Charging Options**: EN 50 desk and wall charger or plug-in transformer.
- ⚡ **Intelligent**. Can be used as mains power supply even while in “charge” mode.
- ⚡ **100 – 240 V Compatible Transformer** incl. 4 plugs for world wide compatibility.

For the optional EN50 table/wall unit please see page 133.



HEINE® LED LoupeLight and Accessories

Ordering Information



LED LoupeLight	
LED LoupeLight for HR/HRP Loupes with i-View loupe mount [01]	C-008.32.246
Detachable Yellow Filter [02]	C-000.32.241
NEW: Polarisation filter P2 for LED LoupeLight and HR Loupes [03] + [04]	C-000.32.535
NEW: Spare polarisation filter P2 for HR Loupes (2 pairs) [04]	C-000.32.536

LED LoupeLight with Power Source for i-View loupe mount (HR / HRP)	
LED LoupeLight with mPack LL and EN50 [05]	C-008.32.261
LED LoupeLight with mPack LL and plug-in transformer [06]	C-008.32.264

HEINE® LED LoupeLight Sets

with plug-in transformer



LED LoupeLight Sets with HR Binocular Loupes			
Magnification and working distance in mm	HR 2.5x/340	HR 2.5x/420	HR 2.5x/520
HR 2.5x optics with i-View on S-FRAME LED LoupeLight Detachable Yellow Filter mPack LL with plug-in transformer Protective lenses 2 sterilisable levers Retaining cord Cleaning fluid Cleaning cloth Carrying case	C-008.32.440	C-008.32.441	C-008.32.442

LED LoupeLight Sets with HRP Binocular Loupes			
Magnification and working distance in mm	HRP 3.5x/420	HRP 4x/340	HRP 6x/340
HRP optics with i-View on S-FRAME LED LoupeLight Detachable Yellow Filter mPack LL with plug-in transformer Protective lenses 2 sterilisable levers Retaining cord Cleaning fluid Cleaning cloth Carrying case	C-008.32.443	C-008.32.444	C-008.32.445